
CRYSTAL WAVE

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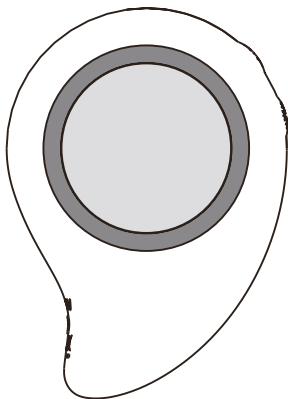
CRYSTAL WAVE

Hot Applicator for Face with Quartz Stone
[CW02P - CW01B] Instructions for Use

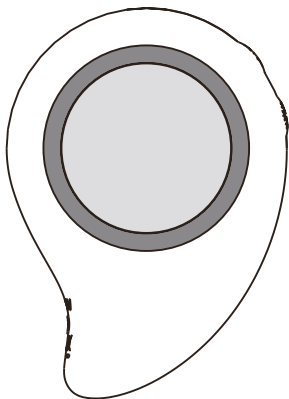Thank you for purchasing **CRYSTAL WAVE** [CW02P - CW01B].
To ensure the safe and proper use of this product, please read the instruction manual carefully and operate according to the guidelines provided.

I. PRODUCT INFORMATION

1. Basic Information

Product Name CRYSTAL WAVE (Hot Applicator for Face with Quartz Stone)

Product Model: CW02P (Amethyst Quartz Applicator)

Product Model: CW01B (Lan Dongling Quartz Applicator)

Dimensions: 87mm * 63mm * 38mm

Charging Parameter: 5.0 V - 1A

Rated Power: 5 W

2. Features

Functions: Heat, vibration, skin cleansing, LED Light Biostimulation

Usage Area: Face

Special Feature: Fully waterproof

Automatic Timer: 15 minutes

Control Mode: Smart button control

3. Scope of Application

3.1 This product can be used on the body.


3.2 Applications: Oily T-zone, enlarged pores, blackheads, acne, dull skin.

4. Features

4.1 Main function: skin cleansing, heat ($42^{\circ}\text{C} \pm 2^{\circ}\text{C}$), vibration massage, LED light biostimulation for skin wellness.

4.2 Two massage modes effectively promote blood circulation, improve the condition of the facial skin, and enhance skin elasticity. The red, green, and blue lights help rejuvenate the skin and improve skin concerns.

II. PRODUCT STRUCTURE

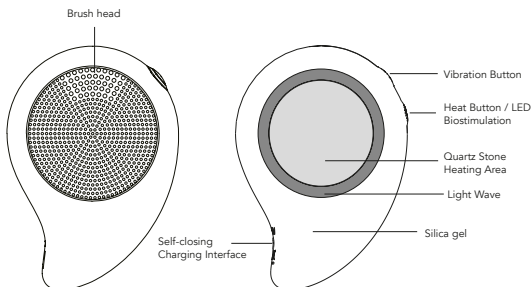
III. HOW TO USE

1. How to Use CRYSTAL WAVE

2. Operating Instructions for CRYSTAL WAVE

2.1 Button Instructions

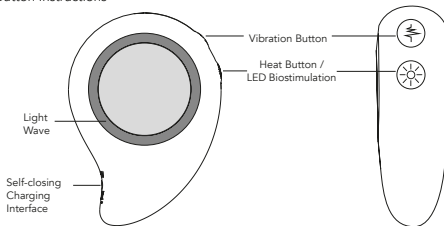2.2 Product Charging (When the product is fully charged, proceed directly to the next step without the need for additional charging.)

2.2.1 Insert the USB cable into the charger socket, then connect it to a household power source. Connect the other end to the device's charging port to start charging. At this stage, the light ring will turn red. Once the battery is fully charged, the red light will turn off.

(Other functions are disabled during charging, and button operation is not available.)

2.3 Operating Instructions

2.3.1

- a) **Power On:** Briefly press the power button to start the device.
- b) **Power Off:** When the device is on, briefly press both buttons (vibration and heat/LED light button)   simultaneously to turn it off.

2.3.2

 Button

 Button for Heat Mode and LED Biostimulation skin treatment

Briefly press the button  ① Green light mode for skin treatment. →

② Blue light mode for skin treatment. → ③ Red light/heat mode. → ④ Off.

2.3.3

 Button

 Button to adjust speed

Briefly press the button  : ① Normal speed → ② Stronger speed

→ ③ Variable frequency speed → ④ Off.

2.3.4 Low Battery: When the battery level is too low, the red light flashes three times, and the device turns off automatically.

2.3.5 Sensor Malfunction: When the sensor is faulty, the red light ring flashes continuously, and the heating function is not operational.

IV. INSTRUCTION MANUAL

1. Vibration for Skin Cleansing

The CRYSTAL WAVE device uses vibration to deeply cleanse the skin, removing impurities, excess sebum, and makeup residues. It offers **three vibration modes**, which can be selected based on skin type and personal preferences:

- **Normal speed:** Recommended for normal skin and daily cleansing.
- **Higher speed:** Suitable for oily skin or deeper cleansing.
- **Variable frequency speed:** Ideal for stimulating circulation and improving skin elasticity.

1.1 USAGE TIME BASED ON SKIN TYPE

- **Normal/Combination Skin:** 1-2 minutes for the entire face.
- **Sensitive Skin:** Use at normal speed for a maximum of 30-60 seconds.
- **Oily Skin:** Up to 3 minutes, with the option to use higher speed on areas prone to impurities.

2. LED Light treatment

The device features three LED lights, each designed for a specific skin treatment purpose:

- **Red Light (Anti-Aging & Regenerating):** Stimulates collagen production, reduces wrinkles, and improves skin elasticity.
- **Green Light (Balancing & Soothing):** Helps reduce redness and dark spots, promoting an even skin tone.
- **Blue Light (Purifying & Antibacterial):** Ideal for acne-prone skin, helping to reduce impurities and inflammation.

2.1 USAGE TIME BASED ON SKIN TYPE

- **Sensitive Skin:** 1-2 minutes per session, up to 3 times a week.
- **Normal Skin:** 3-5 minutes per session, 4-5 times a week.
- **Oily/Acne-Prone Skin:** Up to 5 minutes per session, suitable for daily use.

3. Heat with Red LED Light

The device combines **heated compression** ($42^{\circ}\text{C} \pm 2^{\circ}\text{C}$) with **red LED light**, enhancing the absorption of skin treatments and stimulating microcirculation.

3.1 USAGE TIME BASED ON SKIN TYPE

- **Sensitive Skin:** Maximum **1 minute**, avoiding delicate or irritated areas.
- **Normal/Combination Skin:** **2-3 minutes** across the entire face.
- **Oily Skin:** Up to **5 minutes**, focusing on areas that need more stimulation.

4. Metaphysical Characteristics of Quartz*

The quartz applicator is not just a decoration: each natural crystal is believed to carry energetic vibrations that, according to the principles of crystal therapy, can harmonize energies and bring a touch of well-being and positivity. Ideal for those who want to incorporate spirituality or natural aesthetics into their beauty routine.

- **Amethyst Quartz:** Amethyst is a semi-precious stone that belongs to the quartz family and is used for its relaxing and purifying properties.
- **Lan Dongling Quartz:** Often known as Blue Aventurine, it is considered a stone that helps harmonize emotional balance, bringing tranquility and well-being.


**The information about quartz is for informational purposes only and is in no way intended to diagnose, treat, prevent, or cure medical, psychological, or health conditions. Crystal therapy practice does not replace professional medical consultations, treatments, or care.*

The use of stones and crystals, as indicated in this manual, falls within the scope of complementary and holistic practices aimed solely at personal well-being and energetic balance.

5. Application Techniques and Recommended Movements

To maximize the benefits of the treatment, follow these specific movements while using the device:

- **Neck:** Move from bottom to top to promote skin firmness.
- **Cheeks and Jawline:** Move from the center of the face outward to encourage the drainage of excess fluids toward the lymphatic vessels, achieving a depuffing effect and redefining facial contours.
- **Forehead:** Move from bottom to top to smooth the skin and help prevent wrinkle formation.

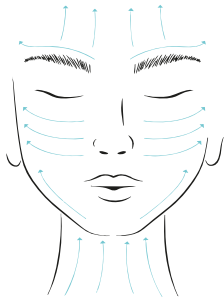By following these recommendations, the treatment will be more effective, leaving the skin looking firmer, more radiant, and relaxed.

V. CLEANING AND MAINTENANCE

1. Cleaning

1.1 This device is waterproof; you can use a cleanser and water to clean it (do not immerse it in water for an extended period).

1.2 Use a damp cloth or a mild detergent to clean the product, then dry it thoroughly.

1.3 Do not use hard brushes, sharp tools, gasoline, or diesel for cleaning.

2. Maintenance

2.1 Store the device in a box and in a safe, dry, and cool place.

2.2 Do not store it in high-temperature areas, near fire, or in direct sunlight to prevent fading or damage.

2.3 Avoid contact with benzene, thinners, or other corrosive liquids.

WARNING: THIS PRODUCT CONTAINS A NON-REMOVABLE LITHIUM BATTERY. DO NOT ATTEMPT TO DISASSEMBLE, PUNCTURE, OR DAMAGE THE BATTERY. AVOID EXPOSURE TO EXTREME HEAT, LIQUIDS, OR OPEN FLAMES. IF THE BATTERY SHOWS SIGNS OF DEFORMATION, LEAKAGE, OR OVERHEATING, DISCONTINUE USE IMMEDIATELY AND CONTACT CUSTOMER SUPPORT. DISPOSE OF THE PRODUCT ONLY THROUGH AUTHORISED CENTRES IN ACCORDANCE WITH LOCAL REGULATIONS.

VI. NOTES

For safety, carefully read the following instructions before use.

1. Use the device according to the instructions in this manual.
2. Do not use parts that are not provided by this manufacturer.
3. Do not use the device while charging.
4. When the battery is fully charged, unplug the charger and do not overcharge the device.
5. Do not store the device in high-temperature areas.
6. Stop using the product in case of damage or other issues and contact customer support.
7. Children must use the device under adult supervision.
8. Individuals with limited physical mobility or reduced sensitivity should use the device under supervision.
9. This product is not intended for medical use.
10. Consult a doctor before use in case of: a) Skin allergies or infections.
11. The device automatically turns off after 15 minutes of use. If further use is needed, restart the device.

VII. TROUBLESHOOTING

Problem	Cause / Solution
DOES NOT TURN ON	<ul style="list-style-type: none"> • Low battery: Fully charge before use. • The product is locked in travel mode: To unlock it, press and hold the power button.
STOPS DURING USE	<ul style="list-style-type: none"> • Low battery: Fully charge before use. • The automatic timer has ended: Restart the device.

Note: If the problem persists, do not use the product and contact customer support (see back cover).

Do not attempt to repair it yourself.

VIII. ENVIRONMENTAL INFORMATION SHEET

The names and content levels of harmful substances in this product are as follows:

Component Name	Sostanza nociva					
	Lead (Pb)	Mercury (Hg)	Cadmium (Cd)	Hexavalent Chromium (Cr(VI))	Polybrominated Biphenyls (PBB)	Polybrominated Diphenyl Ethers (PBDE)
Plastic Components	○	○	○	○	○	○
Leather Case	○	○	○	○	○	○
Electronic Circuit	○	○	○	○	○	○
Power Cable	○	○	○	○	○	○
Hardware Components	○	○	○	○	○	○
Assembled Motor	○	○	○	○	○	○

LEGEND:

This table has been compiled in accordance with the provisions of standard SJ/T 11346.

○: Indicates that the content of the harmful substance in all homogeneous materials of the part is below the required limit, as specified in standards GB4706.1-2005 / GB4706.10-2008.

X: Indicates that at least one content of the harmful substance in the homogeneous materials of the part exceeds the required limit specified in the same standards GB4706.1-2005 / GB4706.10-2008.

IX. QUALITY WARRANTY DECLARATION

Warranty Duration under the Consumer Code

- **Private Customers (B2C)**: In accordance with the Consumer Code, products sold to private consumers are covered by a 24-month legal warranty from the date of delivery. During this period, consumers have the right to request free repair or replacement of the product in case of conformity defects, unless there is evidence of intent or gross negligence on the part of the consumer.
- **Warranty Conditions**: The warranty covers manufacturing defects of the product and the materials used.

The warranty does not cover the following cases:

- **Damage caused by improper use:** The warranty does not cover damage resulting from incorrect use or failure to follow the instructions provided in the manual, including unauthorized commercial use.
- **Tampering or unauthorized modifications:** Any intervention, modification, or repair carried out by unauthorized personnel will void the warranty.
- **Accidental damage:** The warranty does not cover damage resulting from drops, impacts, or any other accidents that may compromise the integrity of the product.
- **Normal wear and tear:** The warranty does not cover the normal deterioration of components subject to wear, such as electrical cables or heating plates, due to continuous use.
- **Damage caused by environmental conditions:** For example, prolonged exposure to heat sources or sunlight.
- **Use with non-compliant electrical power supply:** Damage caused by power surges, lightning, or use with a power source that does not comply with the specified requirements is not covered by the warranty.
- **Cosmetic damage:** Scratches, dents, or other cosmetic damage that does not affect the functionality of the product are not covered by the warranty.
- **Failure to adhere to rest times:** If the product is used beyond the maximum continuous operating time indicated without observing the required rest periods.

To benefit from the warranty, it is necessary to keep the original proof of purchase (receipt, invoice, or e-commerce receipt). In case of defects, the product will be repaired or replaced free of charge at one of our authorized ORGANICS PHARM service centers. To access the warranty service, please contact the following details:

Email: info@organicspharm.com

Phone: (+39) 0445 174 0014

Quality Warranty Card

CUSTOMER NAME		PRODUCT NAME	
CONTACT NUMBER		PRODUCT MODEL	
EMAIL		PURCHASE DATE	

ORGANICS®
P H A R M
VENICE - ITALY

Distributed by: SMÀ srl
Via Bianche,3 - 36010 Carrè (VI) Italy
Hotline: (+39) 0445 1740014 | **Web:** www.organicspharm.com
E-mail: info@organicspharm.com

MADE IN PRC